


This guide only shows adjustments used on a daily basis, please refer to the main end user manual for further information.

Angle Adjustment




To angle Nexus use handset by pressing up and down pads on hand controller. **1** Ensure that the Nexus is free of obstacles which may cause hindrance and damage. Use the angle marker **2** to record the optimum desired position. Do not leave the child unattended. The angle adjustment hand-set controller must be stored safely.

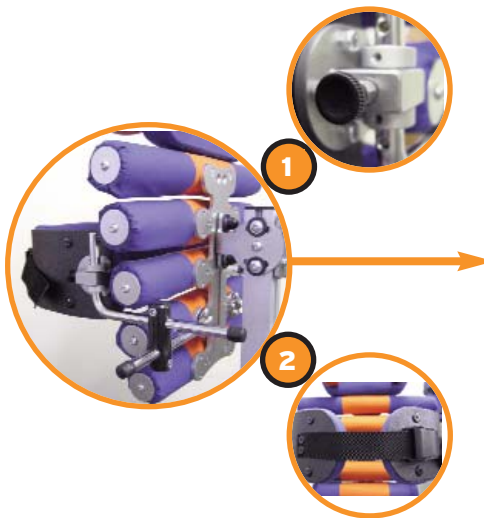

Chest & Pelvic Lateral

To release lateral supports pull plungers **1** and flip support away to the side. Connect buckle across the front of the supports **2**.


- Regularly check chest and pelvic supports are secure and in good condition.

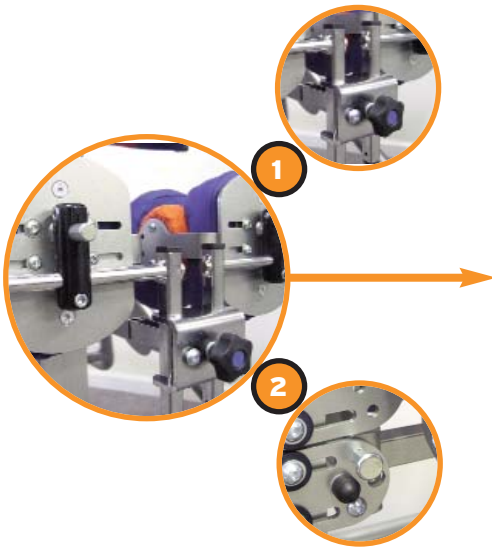

Supine Knee Block

Knee blocks can be removed by undoing the blue knobs by few turns. Like wise to secure the kneeblock slide the knee block to the depth stop and re-tighten the blue knob. ①

Supine knee block fixing stem can be angled up or down when transferring a child by pulling the two plungers out and flipping the stem up or down. ②

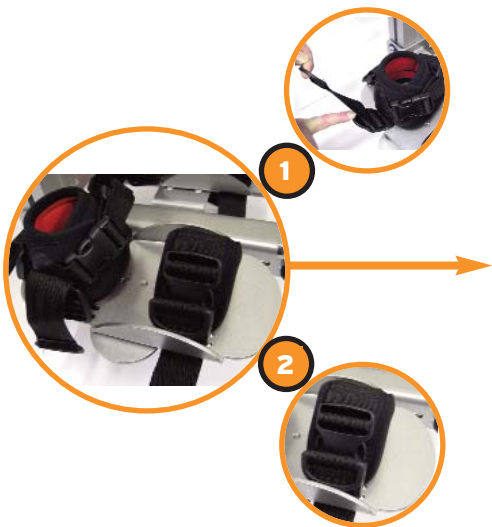

Foot Straps

To fit, loosen webbing at all 3 fixing points, secure cuff around the ankle with buckle fixing and pull each D ring up, to pull webbing tight and to secure the heel. ①

Toe Straps

Fit with buckle fixing adjusting webbing as required. ②

Do **NOT** leave child unattended in the Nexus Stander. Always use all the positioning straps provided.



1



Prone Tray

Loosen the black/orange knobs 1 to achieve the desired tray surface angle adjustment and secure. Fit tray and secure with black wing nuts.



1



Supine Tray

Loosen the orange knobs 1 to achieve the desired tray surface angle adjustment and secure. Fit tray and secure with black wingnuts.



1



Prone Headrest

To fit and remove headrest loosen both black fixing wing nuts 1 on mounting bar.

Important Safeguards

- Do not use this product if you have not received a formal induction from your therapist. As a primary carer please ensure you are aware of your responsibility relating to the safe use of this product
- This product has been prescribed by the Therapist or practitioner who should have adjusted the Nexus to meet the child's need. If this product has not been adjusted by the therapist please do not use it.
- It is imperative that the product log form is completed at the rear of this manual on a regular basis. This form must be used to record setups, changes to the product and scheduled maintenance checks undertaken.
- Please ensure this End User manual is kept in a safe place.
- CARE, MAINTENANCE AND CLEANING procedures must be followed in this manual to ensure client well being.
- It is imperative as a primary carer that you understand the CORRECT USE of the product and also what is INAPPROPRIATE USE.
- Secure the child using the trunk and pelvic laterals, foot straps and knee blocks immediately .
- Always use all of the positioning straps provided.
- Please ensure that the product is used on an indoor even surface.
- It is imperative that children are continuously monitored whilst using this equipment. **NEVER LEAVE A CHILD UNATTENDED IN THIS PRODUCT.**



PO.Box 2125
Wolverhampton, WV6 7FB.
Tel/Fax 0121 523 5004
Tel/Fax 01902 741675.
E-mail: enquire@a2u.co.uk,
Website: www.a2u.co.uk

Serial No : _____

Date : _____

Version 1

